HUMANTECHNIK

lisa
signolux

Wristband receiver

lisa signolux

Wristband receiver

DE	**Bedienungsanleitung**	Seite 4
EN	**User manual**	Page 14
FR	**Mode d'emploi**	Page 24
NL	**Gebruiksaanwijzing**	Pagina 34
FI	**Käyttöohje**	Sivu 44
DA	**Brugervejledning**	Side 54
SV	**Bruksanvisning**	Side 64

REF **A-2646-0**

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The *Lisa-Signolux* »wristband receiver« helps people with hearing impairments compensate for their disability. It allows users to be alerted to events that are detected and transmitted by the *Lisa-Signolux* wireless system.
The »wristband receiver« is part of a medical system and is compatible with all devices of the *Lisa-Signolux* wireless system.

Package contents

Please check to make sure all of the following items are included:

- Wristband receiver
- USB power supply
- User manual
- Charging cable
- 4 country adapters
- Warranty card

If the packaging or the packaging seal (round transparent sticker) is damaged, check the condition of the parts. If any parts are missing or damaged, please contact your retailer or the manufacturer directly.


Immediate safety instructions

- Before commissioning and using the device, read this manual carefully and completely. Observe the other warnings and safety instructions on page 18.
- Keep this user manual in a safe place so that it is accessible to other users at all times. Each main user must be able to access it.
- Keep Lisa-Signolux devices out of the reach of children! This package contains small parts which may present a choking hazard!
- The flashes of light from the receiver devices can trigger epileptic seizures in people with photosensitive epilepsy (PSE).
- After setting up the device, test the signal system to make sure it is working properly!

How It Works / Intended Use

The wristband receiver indicates event signals from the Lisa-Signolux radio transmitters via vibration and colored event symbols.

It also functions as a person call transmitter: When this function is activated, a wireless signal is sent to all receivers within range which will, by default, indicate an alert with the person call icon.

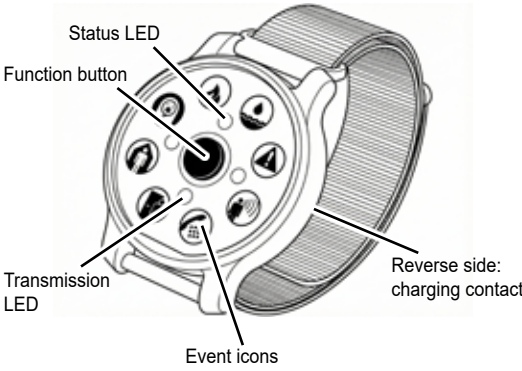Setup

Preparation

If necessary, select the appropriate country-specific adapter for your power outlet and attach it to the USB power supply.

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Battery is charging: The status LED glows red.

Battery is fully charged: The status LED glows green.

When the battery is low, the status LED flashes red. In this case, please connect the receiver to the charger as soon as possible.

Turning on the wristband receiver

When the receiver is off, press and hold the red function button for 2 seconds. All event icons will flash twice, and the receiver will enter standby mode.

Acknowledging / cancelling an active event display:

Press the red function button once briefly to cancel the alert.

Recalling the last received event:

In standby, press the function button once briefly to display the most recently received event again.

Sending a person-call signal:

In standby, press and hold the function button for about 2 seconds – the white transmission LED will flash to indicate transmission. Release the function button as soon as this happens – a person-call signal will be sent to all *Lisa-Signolux* receivers within range.

To switch off the receiver:

In standby, keep holding the function button until all icons light up, then the icons will turn off one by one in sequence.
Release the button before all icons go out. All icons will flash three times and the receiver will switch off.

Automatic pairing

Lisa-Signolux devices can be used together immediately and without prior pairing when delivered.
The displayed icon depends on the type of transmitter. This information can be found in the user manual of the according transmitter.
By default, the person call function is indicated by the person call icon on the receiver units.

Manual pairing

Note 1:
If a transmitter is to be displayed with an event icon of your own choice, this can be done as described in the chapter „*Linking transmitter alerts with event icons*“.
However, this causes the receiver to lose the automatic pairing function and all transmitters previously shown on the receiver must be individually paired.

Note 2:
Should there be repeated false alerts when using the preset transmitter indication, this may hint at a neighbouring Lisa-Signolux installation causing interference. In this case, the problem can be solved by individually pairing all transmitter and receiver units within the setup (see *„Linking transmitter alerts with event icons“*)

Deleting all paired transmitters (re-enabling Automatic pairing):

In standby, briefly press the button first, then immediately press and hold it until all eight icons light up. Release the button - the eight icons will flash once.
This deletes all manually paired *Lisa-Signolux* transmitters from all event icons.

Linking transmitter alerts with event icons

Up to 8 transmitters can be assigned to each event icon. This is achieved by pairing a transmitter with an icon of your choice.
Please keep in mind that a transmitter can only be paired if it is not already stored in the memory of the receiver.
If required, you might first have to delete the transmitter from the memory (see “Deleting all paired transmitters” or „Reset an icon“).

a) In standby, press and hold the function button for more than 5 seconds until all event icons light up in sequence . The icons will cycle as long as the button is held.

b) Let the button go when the desired event icon lights up, then trigger the according *Lisa-Signolux* transmitter.

Test the function after a few seconds by triggering the transmitter again.

Reset an icon

a) In standby, press and hold the function button for more than 5 seconds until all event icons light up in sequence. The icons will cycle as long as the button is held.

b) Let the button go when the desired event icon lights up, then press the button again and hold it until all icons light up. All transmitters paired to that icon are now cleared.

Warnings and safety instructions
This chapter lists further important information for safe handling and other conditions for operation.

⚠ Safety instructions for rechargeable batteries

The wristband receiver has a built-in lithium-ion rechargeable battery. Rechargeable batteries must not be exposed to extremely high or low temperatures or low air pressure during use, storage or transport at high altitudes.
Disposing of the device or the battery it contains by placing it in fire or an oven, or by crushing or cutting it up, may result in fire or explosion.
Exposing batteries to extremely high temperatures may result in an explosion or the release of flammable liquids or gases.
Batteries exposed to extremely low air pressure may explode or release flammable liquids or gases.
Batteries are toxic and must not be swallowed! Keep the device out of the reach of children and pets. In the event of an emergency, seek appropriate medical assistance immediately.

⚠ Information on product safety

Please observe these instructions to avoid damage to the unit and possible invalidation of the warranty.
The device is designed for use in living areas. The unit must not be exposed to direct sunlight for long periods of time and must also be protected from excessive heat, moisture or strong mechanical shocks. Do not drop the unit. Dropping the unit on a hard surface can damage or destroy it.

Opening or tampering with the device will result in the immediate cancellation of the warranty. Use only accessories from Humantechnik.
This device is not suitable for outdoor use. The device must not become damp or wet.

Only use units of the Lisa-Signolux system where electronic devices are permitted. When using in aircraft, it is essential to obtain the prior consent of the flight personnel.

Strong interference signals from third-party devices can impair or completely prevent the function. Do not place potentially interfering devices such as screens and mobile phones closer than 30 cm to Lisa-Signolux devices and, if necessary, their cabling. Important event alerts such as fire warnings may be prevented.

Environmental conditions

Units of the *Lisa-Signolux* system should only be operated in a dry environment within the guideline values for temperature and humidity specified in these operating instructions. Humidity and moisture can lead to failure of the units.

Maintenance and cleaning

The unit is maintenance-free. If it becomes dirty, it can be cleaned with a soft, damp cloth. If necessary, use only commercially available alcohol-based disinfectants. Never use cleaning agents or disinfectants that contain acetone or hydrogen peroxide or that are bleaching, corrosive or abrasive. No liquid or moisture may penetrate the openings of the unit.

Compatible *Lisa-Signolux*-transmitters

The device is compatible with all transmitters of the *Lisa-Signolux* product range:

Universal direct transmitter (A-2655-W)
Universal acoustic transmitter (A-2656-W)
Doorbell-pushbutton (A-2657-W)
Personcall pushbutton (A-2659-W)
Smoke detector Guardion (A-2638-0)
Smoke detector Guardion Interlink (A-2642-0)
Baby transmitter (A-2670-0)
Combined detector CO/Temp (A-2672-0)
Combined detector CO/Temp Interlink (A-2674-0)
Soundmonitor (A-2680-0/W)
Alarmmonitor „Alarmo“ (A-2681-W)
Wristband receiver (A-2646-0)

Specifications

Working range:	10 to 30°C, optimum ambient temperature: 25°C
Operation:	< 90%, non-condensing (transport, storage and operation)
Protection class:	**IP 20**
Radio frequency:	868,425 MHz
Transmission range:	up to 200 m in open spaces, up to 80 m within buildings
Transmission power:	<10 mW
Battery life per charge:	approx. 4 days (standby)
Charging time:	2 hours
Vibration:	Eccentric motor
Display:	all 8 *Lisa-Signolux* event icons
Weight:	30 g (with wristband)
Dimensions:	55 × 43 × 14 mm (receiver/transmitter) 245 mm (strap length)
Battery:	Lithium-ion 190 mAh / 3.6 V

Expected service life of the device:
approx. 10 years (battery not included!)

All serious incidents directly related to this device must be reported to the manufacturer and the relevant local authority!

Technical specifications may be subject to change without prior notice.

Humantechnik's UK authorised representative:

Sarabec Ltd.
15 High Force Road
MIDDLESBROUGH TS2 1RH, United Kingdom

Sarabec Ltd., hereby declares that this device complies with all UK statutory instruments. UK declaration of conformity available from: Sarabec Ltd.

Markings

MD This is the labelling of Humantechnik GmbH for a medical product according to EU 2017/745.

UDI This is the marking for the product identification number to identify the medical device.

SN This is the marking for the serial number to identify the medical device. This is indicated on the device and the packaging.

REF This is the identification for the manufacturer's item number.

This is the marking for the manufacturer of the medical device.

This symbol indicates the production date of the medical device.

CN This is the marking of the country of production of the medical device. This is shown on the device and the packaging.

This symbol indicates distributors for the medical device.

This symbol indicates that the user should read these operating instructions attentively.

This sign indicates warnings that are important for the user as well as important information on use and operational safety.

Temperature range for transport and storage: -5°C to 40°C
Temperature range for operation: 10°C to 30°C, optimum ambient temperature: 25°C

Humidity during transport and storage: < 90% (non-condensing)
Operating humidity: 15% to 90% (non-condensing)

Atmospheric pressure range during transport, storage and operation: 700 hpa to 1060 hpa

NON STERILE This mark indicates that this appliance does not require sterilisation.

CE The CE mark on the device confirms conformity with the European requirements of the Medical Device Regulation (2017 / 745 / EU Medical Device Directive (MDR)).

This device also complies with the following requirements in accordance with European Union directives:

RoHS: Directive 2011/65/EU WEEE: Directive 2012/19/EU
RED: Directive 2014/53/EU

Complete declarations of conformity are available on the Internet at www.humantechnik.com in the service area.

Service, warranty and warranty conditions

The unit has a high level of operational safety. If faults should occur despite proper operation, please contact your specialist dealer or contact Humantechnik GmbH directly.

The warranty covers the repair of the product and returning it to you free of charge. It is recommended that you send in the product in its original packaging, so keep the packaging for the duration of the warranty period.

The warranty does not apply to damage caused by incorrect handling or attempts to repair the unit by people not authorised to do so (destruction of the seal on the unit).

Likewise, damage due to natural forces, such as fire and flooding, is not covered by the warranty. Repairs will only be carried out under warranty if the completed warranty card is returned accompanied by a copy of the dealer's invoice/till receipt.

Always specify the serial number in any event.

Disposal of used electric and electronic units (applicable in the countries of the European Union and other European countries with a separate collection system).

The symbol on the product or the packaging indicates that this product is not to be handled as ordinary household waste but has to be returned to a collecting point for the recycling of electric and electronic units.
You protect the environment and health of your fellow men by the correct disposal of this product.

Environment and health are endangered by a faulty disposal. Material recycling helps to reduce the consumption of raw material.
You will receive further information on the recycling of this product from your local community, your local disposal company or your local dealer.

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